



# belulu® premium GOLD

## User manual (English)

Thank you for purchasing this product.  
Please read this manual carefully before using and keep it well.

### 1 Safety precautions -In order to use safely-

Please read this [Safety precautions] carefully.

For you and others not to hurt. Represents the harm or damage caused by improper handling.

If you give or lend this product to the others, please give this Safety precautions.

#### [Display means]

Display	Means
	Danger
	The wrong way to use, makes people die or risk of serious injury.
	Warning
	The wrong way to use, possible to makes people die or risk of serious injury.
	Attention
	The wrong way to use, damage to the person or property.

#### [Picture display]

	[danger-warning-attention]
	[Prohibited matter]
	[The act of forcing and instruction]

#### Before use

#### ⚠ Danger

- People who use these following medical electronics, do not use this machine
- Implantable medical electronic equipment
- Medical electronic equipment life support
- Medical electronic equipment mounted type
- Other medical electronic equipment
- ※ They may have caused the error.

#### ⚠ Warning

- People who applies to the following
- Acute disease
- Febrile disease
- Tuberculous diseases
- Hemophilia disease
- Abnormal blood pressure
- Malignant tumor
- Infectious disease
- Heart disease
- Angina pectoris
- Multiple sclerosis
- Purpura, hemorrhage
- Skin cancer
- Skin abnormalities
- Pregnancy
- Lactation
- Epilepsy
- Facial neuralgia
- Keloid diathesis
- Atopy
- Excessive tanning
- Purulent inflammation of acne
- Allergic constitution
- Rash experience
- Steroid hormone
- The head and eye have sequel
- Precision eye disease
- Infants and young children
- People who cannot control themselves
- During menstruation
- Other under treatment and abnormality in the body

#### Place-Power

#### ⚠ Attention

- This product is not waterproof. Please do not use or charge in bathroom.
- Do not put it in a place where fire and high temperature & a place where direct sunlight.
- Do not use plug and code with dirt and rubbish.

#### ⚠ Attention

- Please use the specified conditions (voltage, cable etc) Otherwise it will cause fire, electric shock, etc.
- In addition to charging, pull out the cable from the product and outlet. When pulling the plug, do not pull the code.

#### Prepare and use note

#### ⚠ Attention

- Please use on the face.
- Please use the method of the safety precaution.
- Please follow the use of time.

- This product is household electric appliances, do not be used for commercial purposes.
- If the product and the code have any damage, please stop using it.
- Do not touch the machine and the code with wet hands.

#### ⚠ Warning

- If the skin felt abnormal, please stop using it, please consult your doctor or contact our company.
- If the product abnormal, please stop using.

#### ⚠ Attention

- Please follow the following rules when using.
- Please try once, make sure that there is no abnormality after use.
- Take down the metal products.
- Hair removal site, the day can not use.
- ※ May appear skin irritation.
- On the electrical, do not put the long time.
- Please press the all buttons before use.

#### Repair-Maintenance

#### ⚠ Attention

- Do not put it in a place where children can touch.
- Please shut off the power when you clean the product.

- Do not break down or transform the product.


#### ⚠ Warning

- Do not put it in high temperature & direct sunlight
- If you want to throw away the product, please follow the rules of the area where you live.

- Please do not use these things to clean the product, diluent, essential oil, gasoline, containing drugs, cleaning agents, hot water, etc.
- Please do not put the product around the hot machine.
- Please do not fall the product.
- Please do not wash with water.

### 2 Set


※ Accessories and packaging may be have change

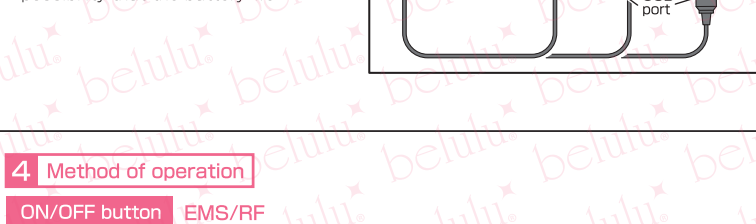

[Product Specification]  
Charging Time: 4~8 hours Input: 5V/0.5A  
Battery type: lithium-ion  
Dimensions: 180x45x35 (mm) Weight: 100g  
[user guide for the supplied AC adapter specification]  
Power Input: 100~240V 50/60Hz 0.2A Power Output: 5V/1A

### 3 CHARGING METHOD


- Please charge battery first.
- When charging, the color of the light will turn red.
- ※ During charging, the machine will generate a little hot, which is normal.
- After charging, the color of the light will turn green.
- In order to protect the battery, even if not used, please charge the battery every month.
- If the operating time is short, there is a possibility that the battery life.

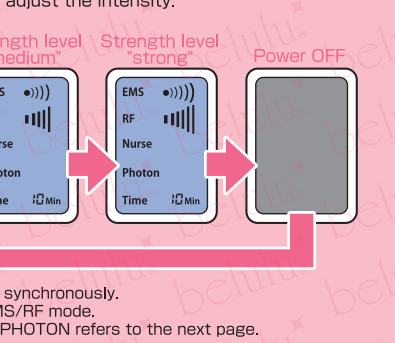

### 4 Method of operation

#### ON/OFF button EMS/RF

- Press ON/OFF button, turn on the power
- The intensity level will be gradually increased by pressing the ON/OFF button each time.
- Turn off the power by press ON/OFF with high intensity.
- Restart the machine to adjust the intensity.

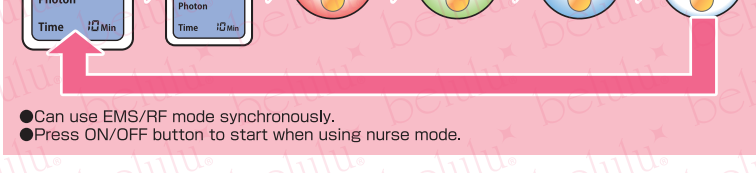

- Can use electroporation and mesotherapy synchronously.
- Press the PHOTON button when using EMS/RF mode.
- Detailed operation about nurse mode and PHOTON refers to the next page.

#### PHOTON button LED & massage

- Press ON/OFF button, turn on the power.
- Turn on the power, and press the PHOTON button to switch and mode: nurse → Red LED → Green LED → Blue LED → Mixed color LED.
- Press the LED button to turn off the LED color under the mixed color LED.
- Press the PHOTON button again to change LED lamp or switch to other models.



- Can use EMS/RF mode synchronously.
- Press ON/OFF button to start when using nurse mode.

#### Combination

EMS/RF hyperthermia + Massage


EMS/RF hyperthermia + LED (PHOTON)

Electroporation and Mesotherapy, will operate at the same time pressing the ON/OFF button.

It is not possible to separately use the Electroporation and Mesotherapy.

#### Time function

Every face care mode is 10min.  
10min later will turn off.

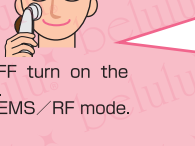


When you press any button, the 10min will restart.

#### [Attention when using]

- Please do allergy test first. ○ Use all of the functions on thin skin, (ex inside the arm)
- If there is an exception, stop using, and consult a doctor.
- When use the LED light, close your eyes, please do not look LED for a long time.
- When use the ion introduction, please use a slide to use beauty liquid. Otherwise dry skin friction is too large. (you can use LED without beauty liquid)
- Keep hand in the ion introduction at the side of the main body of the panel.
- Please straight on the 3D head with respect to the skin. Gently slide slowly when you massage.
- Please take down watches, jewelry and other metal products.
- If you are feel abnormal in use, stop and consult doctor immediately, or contact our company.

You can in accordance with your own skin state to slide the machine.

#### Facial care EMS/RF mode

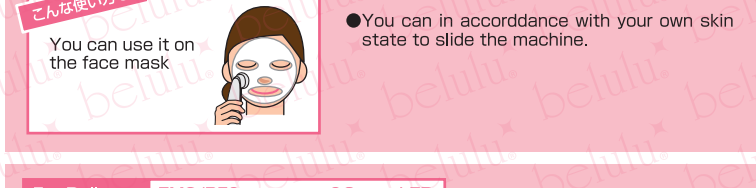


- You can in accordance with your own skin state to slide the machine

#### For dry skin Electroporation/mesotherapy+Nurse



If there is a decrease in beauty liquid, please add, so as to avoid freiction.

#### For Dullness EMS/RF & massage & Green LED

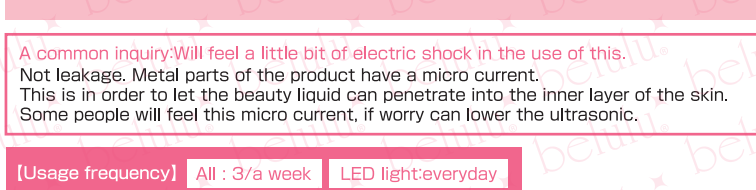


When use the LED light, close your eyes, please do not look LED for a long time.

A common inquiry: Will feel a little bit of electric shock in the use of this. Not leakage. Metal parts of the product have a micro current. This is in order to let the beauty liquid can penetrate into the inner layer of the skin. Some people will feel this micro current, if worry can lower the ultrasonic.


#### [Usage frequency] All : 3/a week LED light: everyday

#### FACE CARE left 2min right 2min Total 4 minutes

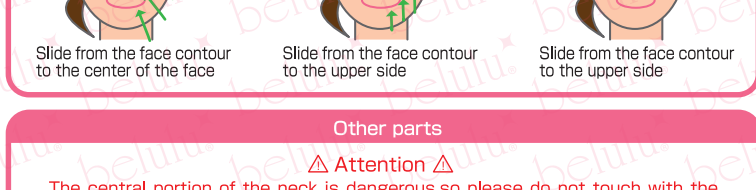

#### Other parts

⚠ Attention ⚠  
The central portion of the neck is dangerous, so please do not touch with the head. Do not put the 3D head into direct contact with the skin above the eyes.

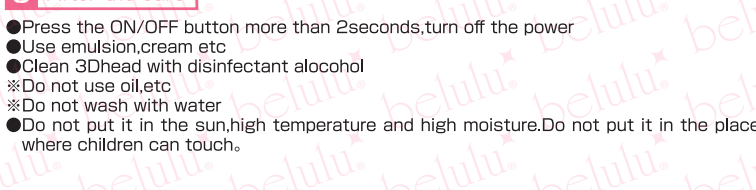

After ordinary slide, slide the same route with circle. Do not put the 3D head into direct contact with the skin above the eyes.

After ordinary slide, slide the same route with circle.

After ordinary slide, slide the same route with circle. The central portion of neck is dangerous, so please do not touch with the 3D head.

### 5 After the care

- Press the ON/OFF button more than 2seconds, turn off the power
- Use emulsion, cream etc
- Clean 3D head with disinfectant alcohol
- ※ Do not use oil, etc
- ※ Do not wash with water
- Do not put it in the sun, high temperature and high moisture. Do not put it in the place where children can touch.

### Q&A

Q How the EMS/RF warm mode run?  
A Under EMS/RF mode, the micro electric current will be transmitted into skin, but won't produce vibration. Here will be short and light stimulation with gentle warm.

Q Could I use it every day?  
A Even though this device brings an excellent cosmetological effect, frequent usage of this device might burden the skin. Therefore, to reach the best result, we recommend using it no more than 3 times a week.

Q Must use beauty liquid ?  
A Toning lotion and cosmetic gel work the same. We strongly recommend our own beauty essence.

Q What is the 3D head made of?  
A On the surface of the 3D head (the part touching the skin), the titanium and 24-karat gold-plated. Not rust, Resistant to scratches, Corrosion-free. You can expect beautiful skin effect.

Q Why the time is always 10min during usage without backward count?  
A The time will be resumed once you press any function button to 10min. So please do not press the other button when using.

Q Weak light (display, LED)  
A Maybe need charge.

If you have any questions Please contact our company

JAPAN  
TEL : 045-227-5121  
E-Mail : shop@kireido.com

CHINA  
TEL : 0771-28752414  
E-Mail : china@kireido.com