

belulu series hikari Manual

Thank you very much for purchasing this product.
Please read this user's manual carefully before using the product.
There is a guarantee certificate at the end of this document,
so please make sure to store it in a safe place.



Warning

- Make sure to use the right voltage, AC adapter, etc. by following the specifications in this document.
- When using this product, do not look directly into the LED light, or point it toward your eyes. Looking directly into the LED may be harmful to the eyes.
- Do not use this product if you find any problems with it, including any damages, faults, deformities, broken cables, etc.
- If you notice any smoke, offensive smell, abnormal heat, etc. coming from the product, turn the power off and stop using it immediately.
- Do not pull on the power adapter cable or touch the plug area with wet hands.
- Do not leave or use this product or the power adapter in places that may be subject to high temperatures, such as stoves, etc.
- Avoid using this product with wet hands, or when it is damp or covered in dirt, etc.
- This product is not waterproof. Never use this product in places with high humidity, such as bathrooms, etc. This may cause burns, electric shocks, or product faults.
- If you notice any abnormal effect on your skin, stop using this product and consult your physician.
- Look at the list below carefully, and do not use this product if any of those conditions apply to you.
- Never attempt to disassemble or reassemble, or do any type of maintenance work on this product by yourself. This may cause the product or its accessories to become damaged or faulty. Any issues caused to the product or its accessories by these actions are not covered by the product guarantee, and the company shall bear no responsibility for them.

- Avoid using this product excessively on sensitive areas of your skin.
- If you experience any pain, etc. during use, move the LED panel immediately away from your skin.
- Do not subject the product to any strong actions such as by dropping or hitting it.
- Do not wash this product with water. If dirty, take a damp cloth, wring it well, and use it to clean the product gently. Do not use solvents such as benzene, thinners, etc. to clean this product.

● Any injuries, faults, damages, etc. caused to the product due to your failure to follow the important notes and warnings stated above, any incorrect use or disassembly/assembly of the product, or any maintenance work done by yourself, are not covered by the product guarantee, and the company shall bear no responsibility for them.

(Danger)

- Acute diseases ● Febrile illnesses ● Blood pressure disorders ● Infectious diseases ● Heart diseases ● Usage of medical electronic equipment (implanted medical devices, such as pacemakers, etc.; electronic life-support equipment, such as heart-lung machines, etc.; or wearable medical electronic devices, such as wearable electrocardiogram devices, etc.) ● History of skin cancer ● Skin disorders caused by diabetes or other diseases ● Currently pregnant ● Currently breast-feeding ● Skin conditions (atopic dermatitis, skin inflammation due to overexposure to the sun, infected and inflamed acne, and sensitive skin due to allergic reactions) ● Dermatitis caused by cosmetics ● Generally not in good physical health ● Usage of long-term doses of steroid hormones drugs, etc.

(Do not apply these parts)

- Malignant tumors ● Cuts ● On bones or joints ● Areas with silicone due to plastic surgery, etc. ● Vocal cords area ● Delicate zones such as hair, around the eyes and ears, mucous membranes, etc.

After using

- Press the power button to turn it off.
- If the product is dirty, use a damp cloth to clean it.
- Avoid storing in places with high temperatures, high humidity, or under direct sunlight, as well as places that are highly dusty or within the reach of children.

This document promises to make the following repairs based on the following conditions.
If you were to encounter a fault within the time period specified below from your date of purchase, you may make a repair request on presentation of this document.

(Repair policy)

1. We would be happy to make repairs if a fault is encountered under normal usage based on the instructions given in this user's manual.
2. Repair requests for faults that have occurred during the product guarantee period may be made by presenting the product and this document.
3. Even if the product is still under the guarantee period, repairs would be chargeable under the following circumstances:
 - a. For any fault or damage caused by the incorrect use, unauthorized repair, or reassembly of the product.
 - b. For any fault caused after purchase during transit, transportation, or due to the product falling to the ground.
 - c. For any fault or damage caused by fire, earthquake, lightning, or any other natural disaster, as well as any public hazard or abnormal voltage.
 - d. If this document is not presented when requesting the repair.
 - e. If this document does not include the customer's name, date of purchase, or if this information has been rewritten over the top.
4. This document is valid in Japan only.

[Warranty]

※Only for the main device only

■No.: ■Name: belulu hikari

■Within 1 year from your date of purchase

■Your name: _____

■Address: _____

■Order #: _____

Beautiful Angel Co., Ltd.

Order date: _____

● When purchasing by mail order, it does not matter if this document does not include the purchase date or store stamp, as customers' information and purchase histories are generally managed by the store. However, your order number will be required when making an inquiry or requesting a repair, so please make sure to make a note of your order number in a safe place.

● This guarantee certificate constitutes our promise to make repairs based on the conditions and the guarantee period specified in this document. Therefore, (for the main device only), this guarantee certificate does not affect customers' statutory rights. Please contact us if you have any queries about repairs to be made outside of the guarantee period.

- Regardless of the number of times that the product has been used, if a repair request is made 3 months after the date of purchase, the customer shall be liable for the return postage fee.
- Outside of the guarantee period, if after diagnosis we consider that repairing the product would be feasible, you may request for such a repair to be executed subject to a repair fee.

We accept a refund only for defective within 7 days.

Please feel free to contact us.

TEL: +81-45-227-5121
E-mail: shop@kireido.com

ZIP: 231-0027

Add: Beautiful Angel Co., Ltd.
Re-tax building 6F Ougicho
2-5-15, Naka-ku
Yokohama, Kanagawa
Japan


belulu hikari Function

420 high luminosity LED beauty lights lead to firm brighter skin!

LED beauty light make skin more active and create a natural beautiful skin.
Its large irradiation area with three-surface, can illuminate whole your face.

The color and function of the LED beauty light.

The basic color is red, blue, and yellow. Each color has different functions.
Choose and care to solve your skin trouble.

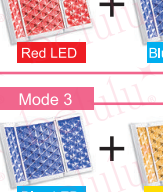


Red LED

Increase collagen, improve wrinkles and flabby skin. It's effective in repairing UV damage skin, and has anti-aging effect.



Blue LED

Remove the active oxygen cause pimples. It can reduce acne and make smooth skin.

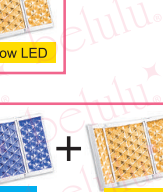


Yellow LED

Promote the process of metabolism, normalize the blood circulation of the skin. Improve the dark spots, acne, pigmentation and other wound problems.


Combination Care

You can irradiate 2 or 3 colors at the same to and save time.

Mode 1

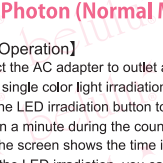


Mode 2

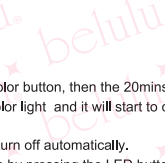

Mode 3



Mode 4

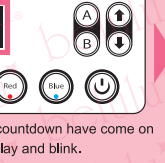


Operation method

LED Photon (Normal Mode)


[Basic Operation]

- Connect the AC adapter to outlet and turn on.
- For the single color light irradiation, press the each color button, then the 20mins countdown have come on the display.
- Press the LED irradiation button to turn on the LED color light and it will start to countdown. (Warning tone will sound once in a minute during the countdown)
- When the screen shows the time is "0", the LED light turn off automatically.
- During the LED irradiation, you can stop the irradiation by pressing the LED button.


[Red LED color light illuminate]

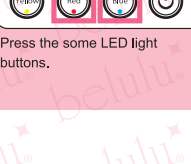

Press the Red LED light button.

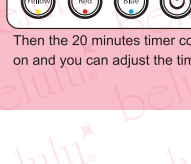


20mins countdown have come on the display and blink.

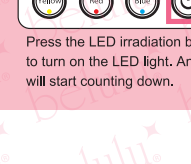


Press the LED irradiation button to turn on the LED light. And it will start counting down.


[Set irradiation time]

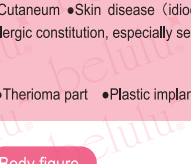

Press the each LED light button and then the 20 minutes timer comes on.

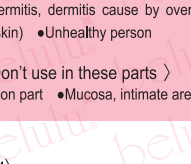

Press the time adjusting button to set the time. You can change within 5-20 minutes.

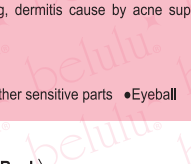

Once set the timer, press the LED irradiation button.

[Simultaneous color lights irradiation]



Press the some LED light buttons.



Then the 20 minutes timer comes on and you can adjust the timer.



Press the LED irradiation button to turn on the LED light. And it will start counting down.

- There may be some changes without notice.

(Don't use followig person)

- Cutaneous ● Skin disease (idiosis dermatitis, dermatitis cause by over tanning, dermatitis cause by acne suppurate, allergic constitution, especially sensitive skin) ● Unhealthy person

(Don't use in these parts)

- Therioma part ● Plastic implanted silicon part ● Mucosa, intimate area and other sensitive parts ● Eyeball


Body figure

(Front)



LED color light irradiation plate

(Back)

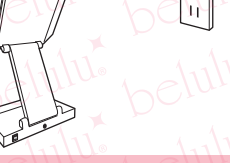

Display and operation keys

Power button

Power port

[Package contents]


- belulu hikari main unit
- Light-cut glasses
- AC adapter

[Preparation for use]

- Insert the power cord to the power port
- Plug in outlet.
- Press the power button behind the device, and it sounds "beep", then turn on.

[Specifications of the product]



- Voltage: AC100-240V(50Hz/60Hz)
- Electricity consumption: 25W
- Weight: 1.4kg
- Size: W:21xL:32xH:7.5cm

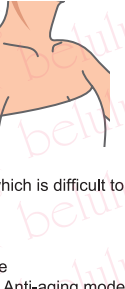

Outlet

- There may be some changes without notice.

belulu hikari also use for these parts

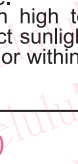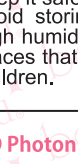

belulu hikari also can be used to above parts which is difficult to care usually.
Please choose proper LED color for you.


For example:

- The ance on the back : Blue LED / Acne mode
- Neck wrinkles : Red LED + Yellow LED light / Anti-aging mode


After using

- After treatment, apply to your skin any beauty serum or cream.
- Turn off the belulu hikari and plug out.
- If the product is dirty, use a damp cloth to clean it. ※Not waterproof
- Please fold the instrument like right figures and keep it safely at proper place.
- Avoid storing in places with high temperatures, high humidity, or under direct sunlight, as well as places that are highly dusty or within the reach of children.

(Figure1) Radiation plate

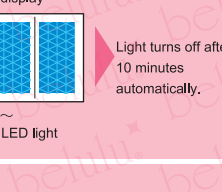


(Figure2)

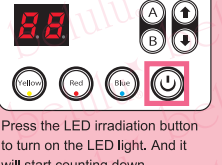


LED Photon (Easy Care Mode)

- Connect the AC adapter to outlet and turn on.
- Choose A or B button, then the countdown have come on the display. You cannot change the time.
- Press the LED irradiation button to turn on the LED color light and it will start to countdown. (Warning tone will sound once in a minute during the countdown)
- When the screen shows the time is "0", the LED light turn off automatically.
- During the LED irradiation, you can stop the irradiation by pressing the LED button.


[Anti-aging mode] (10mins)

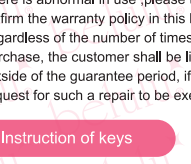

Press the A button.


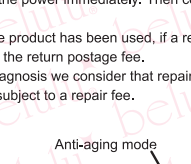


10 minutes countdown comes on the display and blink. You cannot change the time.

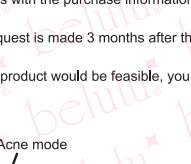


Press the LED irradiation button to turn on the LED light. And it will start counting down.

●Anti-aging care mode procedure

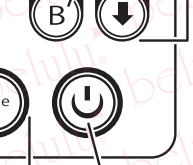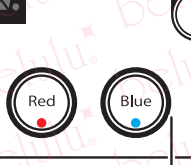


Light turns off after 10 minutes automatically.

[Acne mode] (8mins)

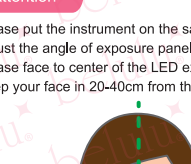

Press the B button.


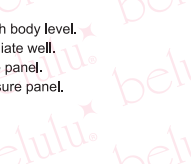


8 minutes countdown comes on the display and blink. You cannot change the time.



Press the LED irradiation button to turn on the LED light. And it will start counting down.

●Acne care mode procedure

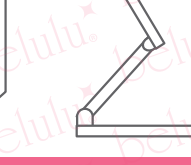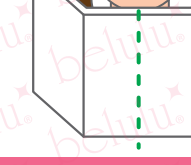

Light turns off after 8 minutes automatically.

After-sale service description

- If there is abnormal in use ,please turn off the power immediately. Then contact us with the purchase information.
- Confirm the warranty policy in this book.
- Regardless of the number of times that the product has been used, if a repair request is made 3 months after the date of purchase, the customer shall be liable for the return postage fee.
- Outside of the guarantee period, if after diagnosis we consider that repairing the product would be feasible, you may request for such a repair to be executed subject to a repair fee.

Instruction of keys



Attention

- Please put the instrument on the same with body level.
- During the LED photons skin care, please wear the light-cut glasses and close your eyes. Looking directly into the LED may be harmful to the eyes.
- Please remove all accessories, watch, jewelry and other precious metal items.
- If you notice any abnormal effect on your skin, stop using this product and consult your physician.

