


# belulu® ~AquaRufa~

## User's manual

Thank you very much for purchasing this product.  
Please read this user's manual carefully before using this product.  
There is a guarantee certificate at the end of this document,  
so please make sure to store it in a safe place.

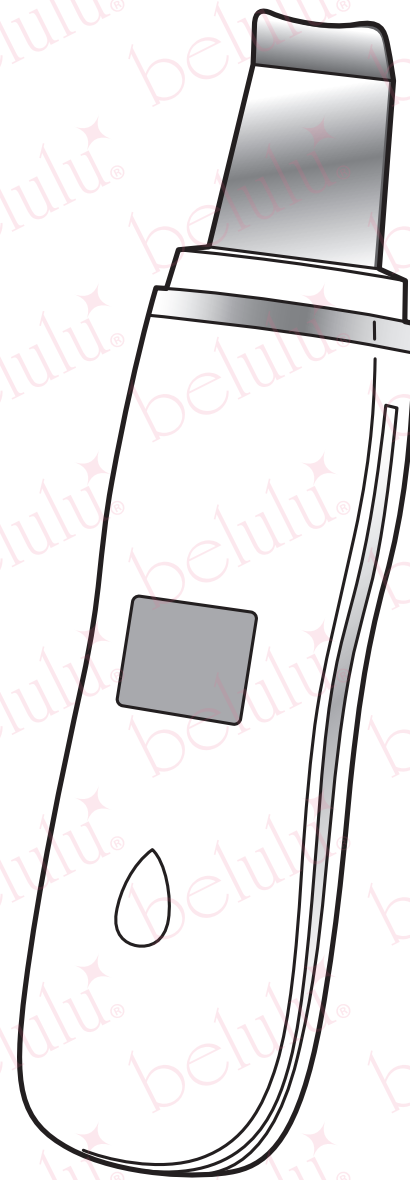

Please check the notice.

### ⚠ Notice

- Make sure to use the correct voltage, AC adapter, etc. by following the specifications in this document.
- When using this product, do not look directly into the LED light, or point it toward your eyes. Looking directly into the LED may be harmful to the eyes.
- Do not use this product if you find any problems with it, including any damages, faults, deformities, broken cables, etc.
- If you notice any smoke, offensive smell, abnormal heat, etc. coming from the product, turn off the power and stop using it immediately.
- Do not pull on the power adapter cable or touch the plug area with wet hands.
- Do not leave or use this product and the power adapter in places that may be subject to high temperatures, such as stoves, etc.
- Avoid using this product with wet hands, or when it is damp or covered in dirt, etc.
- This product is not waterproof. Never use this product in places with high humidity, such as bathrooms, etc. This may cause burns, electric shocks, or product faults.
- If you notice any abnormal effect on your skin, stop using this product and consult your physician.
- Look at the list below carefully, and do not use this product if any of those conditions apply to you.
- Never attempt to disassemble or reassemble or do any type of maintenance work on this product by yourself. This may cause the product or its accessories to become damaged or faulty. Any issues caused to the product or its accessories by these actions are not covered by the product guarantee, and the company shall bear no responsibility for them.

### ⚠ Caution

- Avoid using this product excessively on sensitive areas of your skin.
- If you feel any pain, etc. during use, move the peeling head away from your skin immediately.
- Do not subject the product to any strong actions such as by dropping or hitting it.
- The product is fully charged within approximately 4 to 8 hours. Do not continue to charge it once it is fully charged.
- Disconnect the charger cable from the product when it is not being charged.
- Do not wash this product with water. If dirty, take a damp cloth, wring it well, and use it to clean the product gently. Do not use solvents such as benzene, thinners, etc. to clean this product.
- Any injuries, faults, damages, etc. caused to the product due to your failure to follow the important notes and warnings stated above, any incorrect use or disassembly/reassembly of the product, or any maintenance work done by yourself, are not covered by the product guarantee, and the company shall bear no responsibility for them.

### 〈Do not use people who applies to the following〉


- Acute diseases ● Febrile illnesses ● Blood pressure disorders ● Infectious diseases
- Heart diseases ● Usage of medical electronic equipment (implanted medical devices, such as pacemakers, etc.; electronic life-support equipment, such as heart-lung machines, etc.; or wearable medical electronic devices, such as wearable electrocardiogram devices, etc.)
- History of skin cancer ● Skin disorders caused by diabetes or other diseases ● Currently pregnant ● Currently breast-feeding ● Skin conditions (atopic dermatitis, skin inflammation due to overexposure to the sun, infected and inflamed acne, and sensitive skin due to allergic reactions)
- Dermatitis caused by cosmetics ● Generally not in good physical health ● Usage of long-term doses of steroid hormones drugs, etc.

### 〈Do not use for these parts.〉


- Malignant tumors ● Cuts ● On bones or joints ● Areas with silicone due to plastic surgery, etc. ● Vocal cords area
- Delicate zones such as hair, around the eyes and ears, mucous membranes, etc.

2 Set / Name


※Accessories and packaging may be have changes

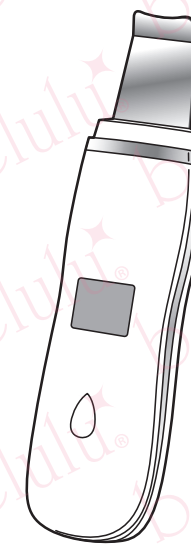

belulu AquaRufa  
main unit

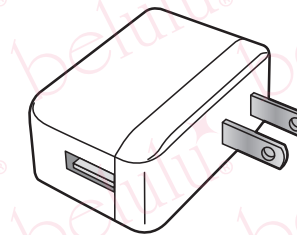

AC Adapter

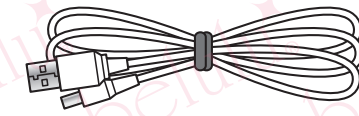

Charge cable

Peeling  
head

Display

operation  
button

Charging  
port

LED light

Peeling  
head

Depending on the setting mode,  
blue or red LED will be lit up.

Ion bar

### 【Main Unit Specification】

Charging: 4~8 hours  
Input: 5V/0.5A  
Battery: Lithium-ion-battery  
Height: 165mm  
Width: 50mm  
Thickness: 20mm  
Weight: 100g

### 【Adapter Specification】

Input: 100~240V 50/60Hz 0.2A  
Output: 5V/1A




### 3 Charging

●The charge notification is lower right in the display. Please charge when the remaining battery is 1 or 0.




\*Depending on the production period, the remaining battery level may be decrease displayed immediately after a full charge, but this is due to a built-in program and there is no problem in use.

- Make sure to recharge the product in the following situations: if you are using the product for the first time after purchase, if you have not used it for a very long time, if its power seems low, or if it seems to work for short periods of time.
- Please turn the power off when you charge.
- Insert the charging cable's connector into the device's charging port and insert the AC adapter into a household power outlet.
- Put a plug into the device, the display will be lit up with beep sounds, and the battery level icon flashing.
- The back light will turn off after a certain time, but it is being charged.
- ※The device may become slightly hotter while charging; however, this is normal.
- The charging time and operating time will vary depending on how frequently the product is used.
- When charging is complete, the battery level icon will stop flashing.
- Then please disconnect the cable.
- This product cannot be used while connected to the charging cable. Make sure to check that the charging cable is disconnected when using the product.
- The most suitable room temperature for charging is 15°C to 30°C. Charging the product at low temperatures may shorten its battery life, so make sure to charge the product at the appropriate temperature. In addition, the product's charging speed may decrease at approximately 0°C or less.
- In order to protect the built-in battery, make sure to charge the product once per month even when you are not using it for long periods of time.
- If the product's operating time seems short even though it has been charged correctly, the cause of this is likely to be the battery life.

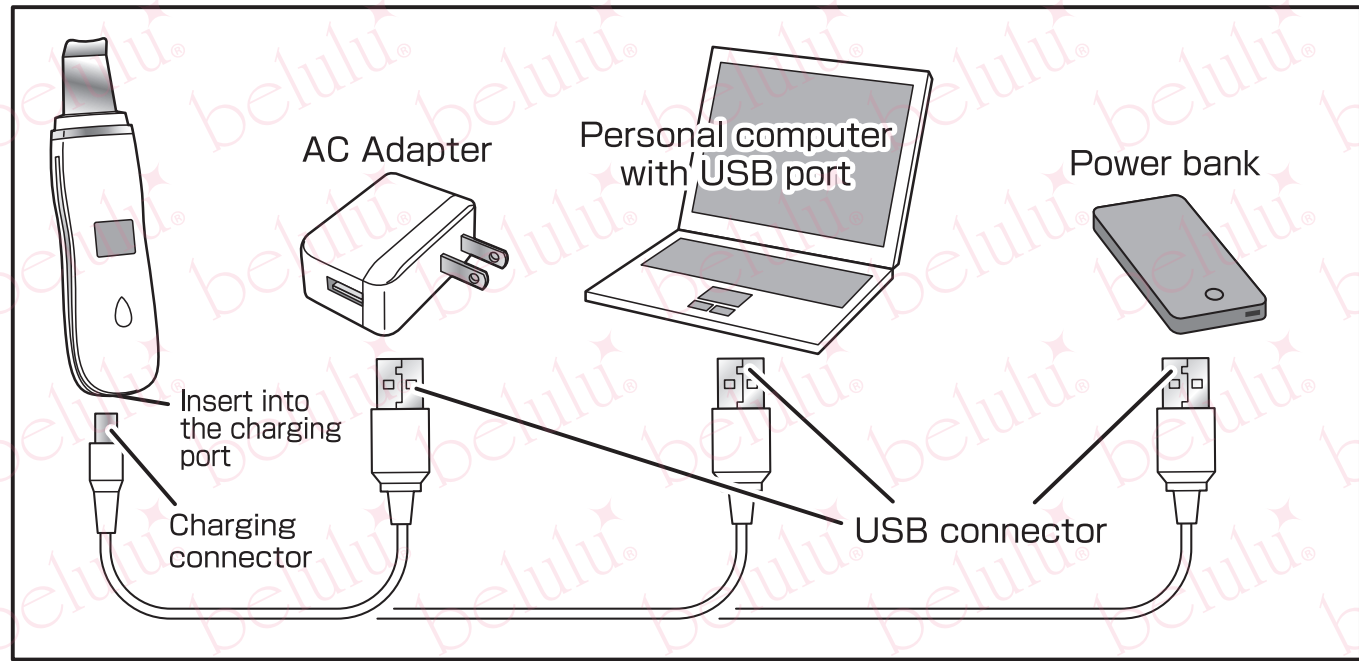

### 4 The display and button operation

Ion mode; + ion mode and -ion mode are displayed respectively

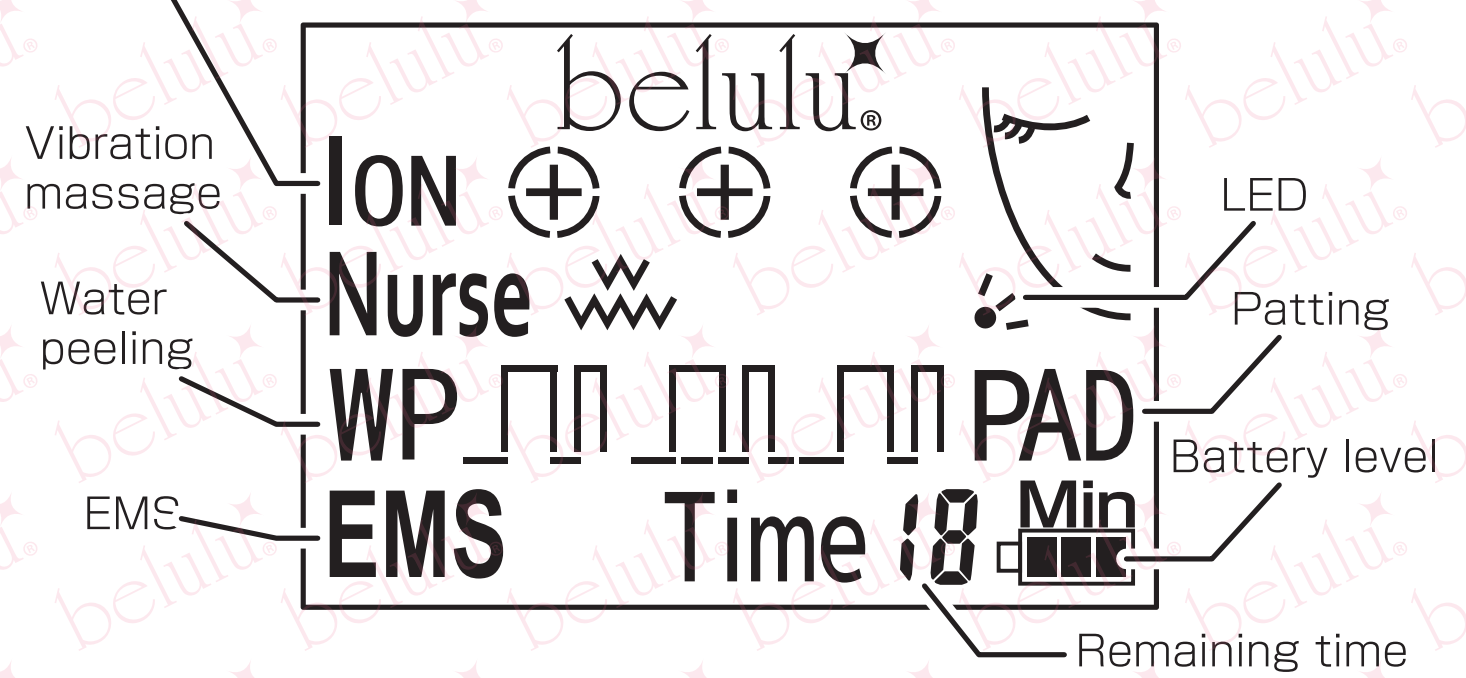


#### Method of operation

##### 【Turn on the power】

- Press the button over two seconds, the display will turn on and stand by state.
- \* For saving power, the display back light will turn off after a certain time. Lights up again when the button is operated.

##### 【Standby state】

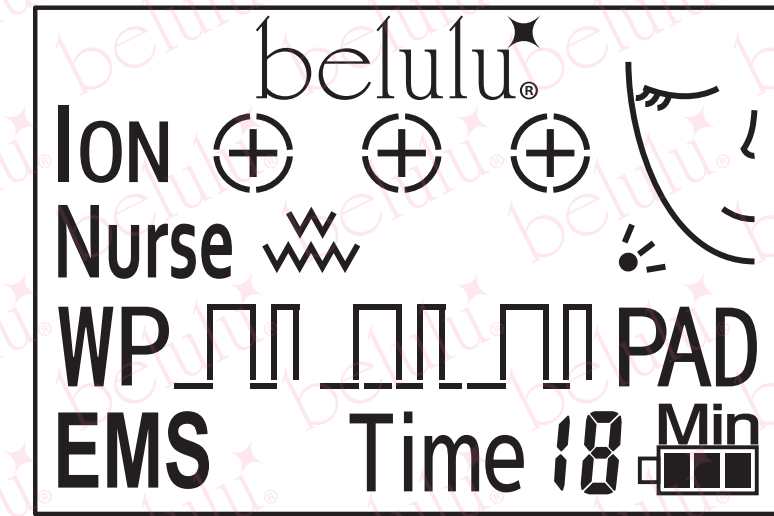

This product operates several functions at same time for each mode. Each function cannot be used separately.


### 【Cleansing mode】

- Press the button once, the cleansing mode will be operated.

#### 【Cleansing mode operation screen】




#### 【Functions that operate at the same time】

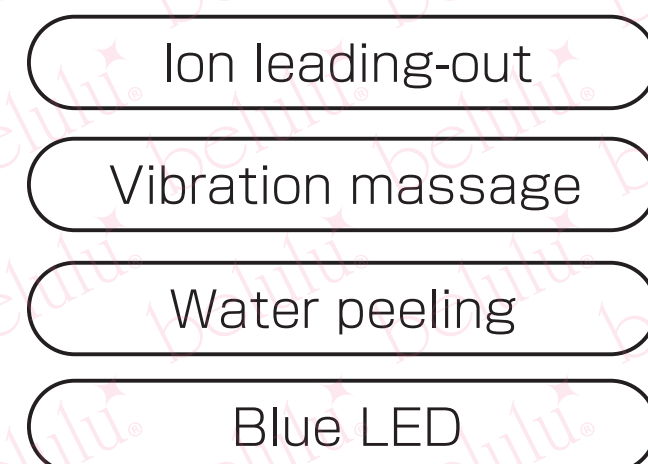


### 【Lift-up mode】

- Press the button twice, the cleansing mode will be operated.  
If the cleansing mode is already operated, the lift-up mode will be starting by pressing the button once.

#### 【Lift-up mode operation screen】

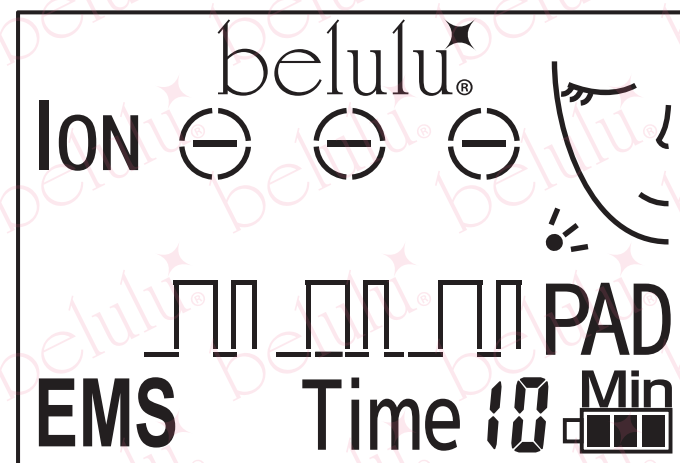

#### 【Functions that operate at the same time】

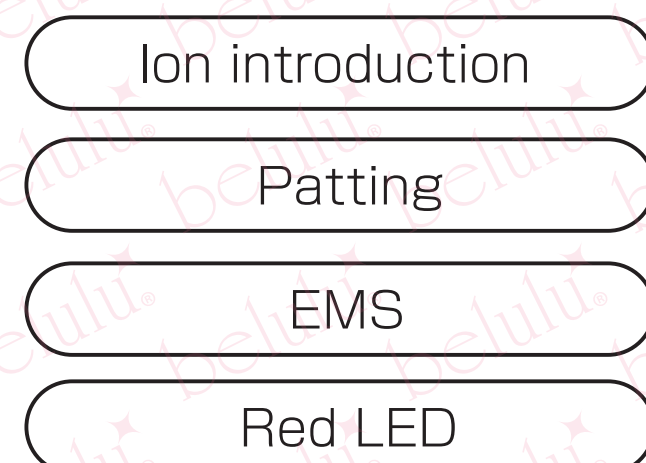

## 5 Before you use

Make sure to do an operation check and patch test before you use.

### Operation check


#### 【Water peeling mode】

1. Refer to the 【Button Operation】 page and set the cleansing mode.
2. Wet your skin with water and gently slide the peeling head against the skin. If a fine mist comes out from the tip of the peeling head, it is operating normally.

It can also be confirmed by dropping water directly on the peeling head.

- \* The mist may be difficult to see due to room brightness and background color.

- \* Pay attention to the direction of the peeling head.  
(Please refer the picture below.)

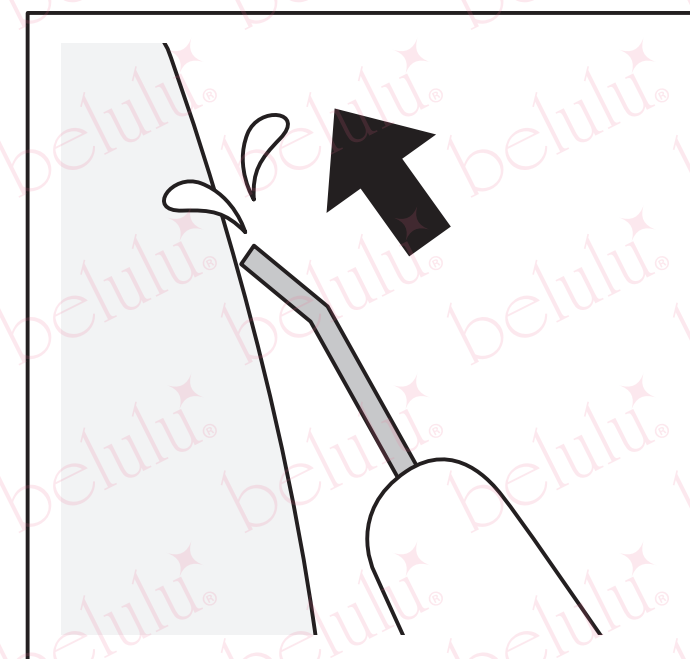

### patch test

- Make sure to do an operation check and patch test before you use.
- Wet the inside of the arm with water and try cleansing mode.
- Please wipe off the water inside the arm, apply beauty gel or your gel, and try the lift-up mode.
- Once you have verified that you have had no adverse reactions, you can move onto your face care. If you do experience any abnormality on your skin, do not use the product, and make sure to consult our company or your physician.



## 6 Function and effect

### 【Water peeling】

The device removes and cleans the moisture on the skin by sonic vibration while blowing away dirt, sebum, old keratin, the back of the pores etc. There is no pain or heat, and it is gentler to your skin than scrubs and pore packs!

### 【【Ion leading-out】】

For expelling pore dirt, grease off to deep clean skin, make the skin to rejuvenate.

### 【Vibration massage】

Vibration can make the nutrition and moist penetrate to the deep skin for effective absorption and moist replenishment, accelerate skincare absorption. Also help to improve sluggish blood circulation.

※Ion introduction works at the same time. You can't use separately.

### 【Blue LED】

Bacteria in pimples contain Porphyrins; able to kill Porphyrins due to the high affinity. The blue light furthermore has calming action which is very effective for hypersensitivity.

### 【【Ion introduction】】

Can make the nutrition and moist penetrate into the deep skin for effective absorption and moist replenishment, accelerate skincare absorption

※Vibration massage works at the same time. You can't use separately.

### 【Patting】

Also gives your skin firmness with a micro current and tapping effect. Encourage turnover and smooth skin.

And it may feel nice lift-up effect!

### 【EMS】


The EMS stimulate your facial muscle and makes your face look smaller. Also activating metabolism of skin, and collagen growth.

### 【Red LED】

Will stimulate the production of Collagen is an essential protein used to repair damaged tissue and to replace old tissue. Best for fine lines and will also help to reduce large pores.

## 7 Let's start!

### Basic care

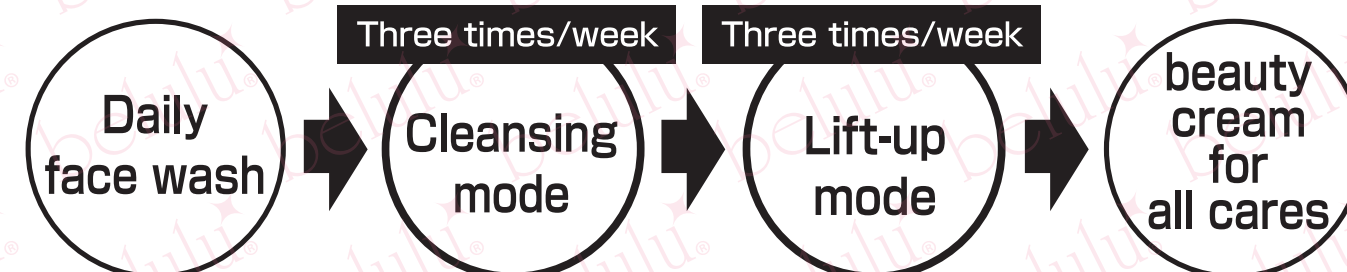

### Caution

- Wash your face as usual and wipe away any moisture. When using water peeling it is required to wet your face, but we recommend wiping off the water once after you wash your face.
- Do not use on dry skin. It may cause skin trouble.
- Please keep moving the peeling head and do not keep on the same part of the skin.
- If the peeling head is difficult to slid, add water, serum, beauty gel, and lotion each time.
- Do not press the peeling head strongly against your skin.
- Remove precious metals such as watches and jewelry when using.
- If you feel pain or other abnormalities during use, stop using it immediately and consult a doctor.

### 【What to prepare】

- Water or cleansing lotion
  - …For Ion leading-out & Water peeling
- Beauty serum, gel or facial lotion
  - …For Ion introduction. We recommend thick beauty serum which contains active ingredients blended, especially 'belulu b2 bottle'.
- Milky lotion or beauty cream
  - …For all cares

## Cleansing mode

10 minutes / day, 3 times / week.

Remove your makeup and wipe off moisture.

Remove pore dirt, sebum, and old keratin by fine mist and ion-leading out.



- ① Wet the entire face with water or wiping lotion.  
Do not use on dry skin.

If the skin is dry, the skin will be damaged by friction, so please wet it with plenty of water.



- ② Refer to 【Button operation】 to activate cleansing mode.  
Grip the micro current induction sidebar, gently touch the tip of the peeling head against your skin and slide it across the face push the peeling head softly.

\*Be careful do not to look at the LED light for a long time.  
If you feel dazzling, close your eyes.

## How to move the peeling head -Cleansing-

This is the basic way to move.  
Such as focusing on the areas that you care,  
but please arrange according to your skin condition.

### Care around the cheek

- ① Starting from the face line, move it slowly to the cheekbones.  
② Move slowly from the center of the face to the outside, shifting little by little thin the range of the cheekbones.



### Care around the T zone

- ① From bottom to top to follow the shape of the nose. The nose circumference cares tilting the nose left and right.  
② Move the forehead from the top of the eyebrows to the hairline.



### Care around the mouth

- ① Slowly move from the center to outside.  
② Pay attention to the corners of the mouth where dullness and dirt tend to accumulate. It is effective to move it little by little in various directions.





## Lift-up mode

10 minutes / day, 3 times / week.

Wipe off water remaining on the face after cleansing mode.

Moisturizing and lift-up with micro current and ion introduction.



- ① Please apply beauty serum, gel, lotion, etc. to your face. (If the gel has viscosity, it can be put on the back of the peeling head.)

\*Do not use on dry skin.



- ② Activate the lift-up mode by referring to 【Button operation】. Grip the micro current induction side bar, gently touch the back of the peeling head against your skin and slide it across the face to pull the peeling head.

## Q&A

### Can I use my own lotion and serums?

You can use your favorite lotion, essence or beauty gel when patting. However, we recommend that you use a gel-like one that does not contain any additives and stays on the skin, or one that is thick. Face pack is also recommended.

## How to move the peeling head -Lift up-

This is the basic way to move.

Such as focusing on the areas that you care, but please arrange according to your skin condition.

### Cheek

- ① Starting from the face line, move it slowly to the cheekbone in a few times.
- ② Move slowly the peeling head from the center of the face to the outside, shifting little by little within the range of the cheekbones.



### Forehead, Nose

- ① Move the peeling head from the side of the nose to the top with the image of pulling up the slack line and the slack of the cheeks.
- ② Move the forehead from the top of the eyebrows to the hairline. Move to pull up both ends of the forehead from around the corner of the eyes.



### mouth, chin

- ① Slowly move from the center to outside.
- ② Pay attention to the corners of the mouth where slackening. It is effective to move it little by little in various directions.



## 8 After care

- Turn off the power.
- Blend in the rest of the gel, facial lotion, or beauty serum on your face.
- Seal it with a milky lotion or beauty cream to compensate for the fat content.
- After use, wiping off the gel, lotion, and serum remaining on the peeling head for clean. \* Only the peeling head is water resistant. The parts other than the peeling head are not completely water resistant, so do not wash them with water.
  - ※ Do not use solvents such as benzine, thinners, etc. to clean this product.
- Avoid storing in places with high temperatures, high humidity, or under direct sunlight, as well as places that are highly dusty or within the reach of children.

## 9 Q&A

<p>The device doesn't work.</p>	<ul style="list-style-type: none"><li>● <b>It may be out of charge.</b> The rechargeable battery has the peculiar attribute of losing power a little at a time even when not turned on, and there may be times when it has no charge. Recharge for one hour and try again.</li><li>● <b>Press the ON/OFF button over two seconds, the display will turn on and standby state.</b></li><li>● <b>After a while, press the operation button again.</b> Operation buttons may be temporarily disabled immediately after the power off, but this is not a malfunction. Wait for about three seconds and then press the operation button again.</li></ul> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>
<p>Not charged, not stable charging.</p>	<ul style="list-style-type: none"><li>● <b>The cable is connected securely?</b> It can't be charged if it isn't energized. Please refer the 'Charging' and check the charging port and connect securely.</li></ul> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product, cable or battery into consideration, so please contact us.</p>
<p>Ion leading-out and ion introduction don't work.</p>	<ul style="list-style-type: none"><li>● <b>Check the operation.</b> Due to Ion leading-out and Ion introduction's feature of sending a weak electric current into the skin, the oscillations or stimulus cannot be heard.<ul style="list-style-type: none"><li>※ You can hear the Vibration massage sounds. Please refer to '5. Before you use' and 'Ion leading-out', 'Ion introduction' and check the operation.</li></ul></li></ul> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>



<p><b>Water peeling mist cannot be seen</b></p>	<ul style="list-style-type: none"> <li>●<b>Increase the amount of water.</b> If the skin is dry, the mist will not rise. Please use plenty of water. The water peeling mist is so fine that it may be difficult to see depending on the ambient brightness and background color.</li> <li>●<b>Check the mode.</b> Please refer to <b>【Button Operation】</b> and check if it is in cleansing mode.</li> <li>●<b>Check the remaining battery.</b> If the remaining battery is low, the mist may not increase due to insufficient output. If the remaining battery is one or zero, please charge the device.</li> </ul> <p><b>【Still have problem?】</b> We have taken early failures and breakdowns of this product into consideration, so please contact us.</p>
<p><b>When in use, my hands will tingle as if it is short circuiting.</b></p>	<ul style="list-style-type: none"> <li>●<b>This is not a short circuit.</b> A tiny amount of electricity is flowing through the metal handle of this device, and it is meant to pass through your skin's barrier to permeate its beautifying elements deep in the skin, so it is possible depending on the person that this tiny electric current will be felt. This small electric current is harmless to your body, but if you are concerned, use it by turning down the strength level of the ultrasonic waves.</li> </ul> <p><b>【Still have problem?】</b> Please contact us. Within 7 days from arrival.</p>
<p><b>During use, there is a "beep" and it powers off.</b></p>	<p><b>It powers off after 10 minutes automatically.</b></p>

**Q. What is the material of the peeling head part?**

**A.** It is stainless.  
It is strong against cracks, corrosion, or rusting, and it is hard for it to trigger allergies. It may not be difficult for it to trigger allergies for everyone. Please do a patch test.

**Q. My skin is stiff, so right after I wash my face, I want to use Ion leading-out after applying the facial lotion.**

**A.** Ion leading-out has the function of absorbing the dirt not removed from washing your face, but if you apply facial lotion right after washing your face, the dirt will remain and clog your pores, so there will be no effect from the Ion leading-out and Ion introduction after that. Stiff skin is caused by chronic dryness of the inner part of the skin, so you will need to thoroughly moisturize it. It will hurt at first, but strictly adhere to carefully removing the dirt by Ion leading-out and then moisturize by Ion introduction. Continual use moisturizes the inner part of the skin, improving stiff or dry skin.

**Q. Is it safe to use if I'm pregnant?**

**A.** You can use it even if you are pregnant, but your skin is in a delicate condition from the hormonal imbalance during pregnancy, so depending on the timing when it is used, it is possible in rare cases that it will harm your skin. It is recommended that you make sure to consult your physician and use it while observing how it affects your skin.  
(The beauty device's functions or the ingredients from the beauty serum have no effect on a mother's body or a fetus, so please don't worry about that.)

**Q. Even after using it, I can't feel any of its effects.**

**A.** How the effects materialize is different depending on the person. This device was manufactured with the intention of it being used at home, so it does not have the same power as from a beauty salon, and there is no immediate effect. By continuing to use it every day, you will start to see its effects.

**Q. Should I take time every day to perform this care routine?**

**A.** This device's built in functions demonstrate highly beautifying effects, but overuse will put a burden on your skin. So as not to burden your skin during use, abide by the frequency of use as listed in this manual. Frequent use every day will not produce immediate effects, it is important to carefully continue to use it on a daily basis.

If you have any questions or concerns,  
please feel free to contact us.

**JAPAN:** [info@kireido.com](mailto:info@kireido.com)

**CHINA:** [china@kireido.com](mailto:china@kireido.com)